

REFERENCES

1. Schall, T.J. et al. (1988) J. Immunol. 141:1018.
2. Schall, T.J. (1991) Cytokine 3:165.
3. Schall, T.J. et al. (1992) Eur. J. Immunol. 22:1477.
4. Heeger, P. et al. (1992) Kidney Int. 41:220.
5. Schall, T.J. et al. (1990) Nature 347:669.
6. Taub, D.D. and J.J. Oppenheim (1993) Cytokine 5:175.
7. Loetscher, P. et al. (1996) J. Immunol. 156:322.
8. Bischoff, S.C. et al. (1993) Eur. J. Immunol. 23:761.
9. Kameyoshi, Y. et al. (1992) J. Exp. Med. 176:587.
10. Lin, C.L. et al. (1998) Eur. J. Immunol. 28:4114.
11. Mattoli, S. et al. (1995) Biochem. Biophys. Res. Comm. 209:316.
12. Fine, J.S. et al. (2001) Inflammation 25:61.
13. Cross, A.K. and M.N. Woodroffe (1999) J. Neurosci. Res. 55:17.
14. Turner, L. et al. (1995) J. Immunol. 155:2437.
15. Bacon, K.B. et al. (1995) Science 269:1727.
16. Appay, V. et al. (2000) Int. Immunol. 12:1173.
17. Appay, V. et al. (1999) J. Biol. Chem. 274:27505.
18. Fischer, F.R. et al. (2001) J. Immunol. 167:1673.
19. Zhang, Y. et al. (2003) FASEB J. 17:734.
20. Proudfoot, A.E.I. et al. (2003) Proc. Natl. Acad. Sci. USA 100:1885.

Human CCL5 Immunoassay

Catalog Number: SEKH-0061

For the quantitative determination of human CCL5 concentrations in cell culture supernates, serum, and plasma.

For research use only. Not for use in diagnostic procedures.

Country | Company: China | Beijing Solarbio Science & Technology Co.,Ltd
Address:NO.85A, Liandong U Valley, Tongzhou District, Beijing, P.R.China.
Tel: 86-10-56371241 Fax: 86-10-56371282 E-mail: service@solarbio.com

TABLE OF CONTENTS

SECTION	PAGE
BACKGROUND.....	01
PRINCIPLE OF THE ASSAY.....	01
TECHNICAL HINTS AND LIMITATIONS.....	02
PRECAUTIONS.....	02
KIT COMPONENTS& STORAGE CONDITIONS.....	03
OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED.....	04
SPECIMEN COLLECTION & STORAGE.....	04
REAGENTS PREPARATION.....	04
ASSAY PROCEDURE.....	06
CALCULATION OF RESULTS.....	06
PERFORMANCE CHARACTERISTICS.....	08
REFERENCES.....	10

LINEARITY: To assess the linearity of the assay, three samples were spiked with high concentrations of CCL5 in various matrices and diluted with the appropriate Sample Diluent to produce samples with values within the dynamic range of the assay.

Dilution ratio	Recovery(%)	Citrate plasma	Cell culture supernatants
1:2	Average% of Expected	98	104
	Range(%)	89-107	96-113
1:4	Average% of Expected	102	106
	Range(%)	94-114	99-115

Performance Characteristics

SENSITIVITY: The minimum detectable dose was 3.9pg/mL.

SPECIFICITY: This assay recognizes both natural and recombinant human CCL5. The factors listed below were prepared at 100ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

Factors assayed for cross-reactivity

Recombinant human	Recombinant mouse	Recombinant porcine
ANG	GM-CSF	PDGF
CNTF	IL-1 α	TGF- β 1
β -ECGF	IL-1 β	TGF- β 2
EGF	IL-3	
FGF acidic	IL-4	
G-CSF	IL-5	
GM-CSF	MIP-1 α	
IL-1 α	MIP-1 β	
TGF- α	SCF	
TNF- α	TNF- α	

REPEATABILITY: The coefficient of variation of both intra-assay and inter-assay were less than 10%.

RECOVERY: The recovery of CCL5 spiked to three different levels in four samples throughout the range of the assay in various matrices was evaluated.

Recovery of CCL5 in two matrices

Sample Type	Average % of Expected Range(%)	Range(%)
Citrate plasma	93	86-102
Cell culture supernatants	97	89-106

BACKGROUND

RANTES (Regulated upon Activation, Normal T cell Expressed and presumably Secreted), also known as CCL5, is a member of the "CC" subfamily of chemokines. It plays a primary role in the inflammatory immune response via its ability to chemoattract leukocytes and modulate their function. The cDNA for RANTES was initially discovered by subtractive hybridization as a T cell specific sequence. Human RANTES cDNA encodes a highly basic 91 amino acid (aa) residue precursor polypeptide with a 23 aa hydrophobic signal peptide that is cleaved to generate the 68 aa mature protein. Human RANTES exhibits approximately 85% homology with mouse RANTES at the deduced aa level.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for CCL5 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any CCL5 present is captured by the coated antibody after incubation. Following extensive washing, a biotin-conjugate antibody specific for CCL5 is added to detect the captured CCL5 protein in sample. For signal development, horseradish peroxidase (HRP)-conjugated Streptavidin is added, followed by tetramethyl-benzidine (TMB) reagent. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Solution containing sulfuric acid is used to stop color development and the color intensity which is proportional to the quantity of bound protein is measurable at 450nm.

DESCRIPTION



TECHNICAL HINTS AND LIMITATIONS

1. This Solarbio ELISA should not be used beyond the expiration data on the kit label.
2. To avoid cross-contamination, use a fresh reagent reservoir and pipette tips for each step.
3. To ensure accurate results, some details, such as technique, plasticware and water sources should be emphasized.
4. A thorough and consistent wash technique is essential for proper assay performance.
5. A standard curve should be generated for each set of samples assayed.
6. It is recommended that all standards and samples be assayed in duplicate.
7. Avoid microbial contamination of reagents and buffers. Buffers containing protein should be made under aseptic conditions and be prepared fresh daily.
8. In order to ensure the accuracy of the results, the standard curve should be made every time.

PRECAUTIONS


The Stop Solution suggested for use with this kit is an acid solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

DESCRIPTION

regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Typical data using the CCL5 ELISA

Standard (pg/ml)	OD.	OD.	Average	Corrected
0	0.048	0.051	0.050	-
7.81	0.190	0.188	0.189	0.139
15.6	0.276	0.273	0.275	0.225
31	0.411	0.406	0.408	0.359
62.5	0.652	0.645	0.648	0.599
125	1.059	1.047	1.053	1.003
250	1.725	1.705	1.715	1.665
500	2.810	2.777	2.794	2.744

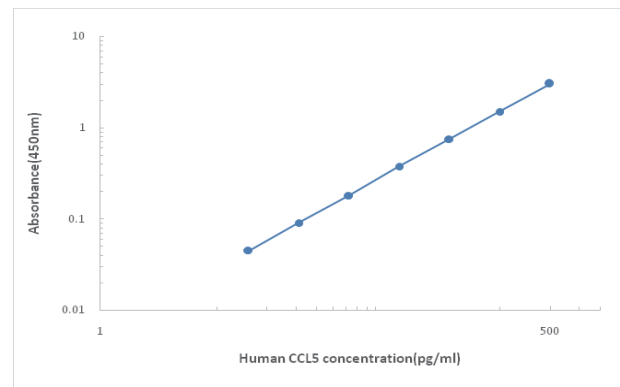

Representative standard curve for CCL5 ELISA.

ASSAY PROCEDURE

Prepare all reagents and standards as directed.

Wash the plate 3 times before assay.



Add 100µl standard or samples to each well, incubate 90 minutes, 37°C.



Aspirate and wash 4

Add 100µl working solution of Biotin-Conjugate anti-human CCL5 antibody to each well, incubate 60 minutes, 37°C.



Aspirate and wash 4

Add 100µl working solution of Streptavidin-HRP to each well, incubate 30 minutes, 37°C.



Aspirate and wash 5

Add 100µl Substrate solution to each well, incubate 15 minutes, 37°C.
Protect from light.



Add 50µl Stop solution to each well. Read at 450nm within 30 minutes.

CALCULATION OF RESULTS

1. The standard curve is used to determine the amount of specimens.
2. First, average the duplicate readings for each standard, control, and sample. All O.D. values are subtracted by the mean value of blank control before result interpretation.
3. Construct a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
4. The data may be linearized by plotting the log of the CCL5 concentrations versus the log of the O.D. and the best fit line can be determined by

KIT COMPONENTS & STORAGE CONDITIONS

PART	SIZE	STORAGE OF OPENED/ RECONSTITUTED MATERIAL
Microwell Plate - antibody coated 96-well Microplate (8 wells x12 strips)	1 plate	Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2 – 8°C**
Standard - lyophilized, 2000 pg/ml upon reconstitution	2 vials	Store at 2-8°C **for six months
Concentrated Biotin-Conjugated antibody(100X) - 120 ul/vial	1 vial	Store at 2-8°C ***for six months
Concentrated Streptavidin-HRP solution(40X) - 300 ul/vial	1 vial	Store at 2-8°C ***for six months
Standard /sample Diluent - 16 ml/vial	1 bottle	Store at 2-8°C ***for six months
Biotin-Conjugate antibody Diluent - 16 ml/vial	1 bottle	Store at 2-8°C ***for six months
Streptavidin-HRP Diluent - 16 ml/vial	1 bottle	Store at 2-8°C ***for six months
Wash Buffer Concentrate (20x) - 30 ml/vial	1 bottle	Store at 2-8°C ***for six months
Substrate Solution - 12 ml/vial	1 bottle	Store at 2-8°C ***for six months
Stop Solution - 12 ml/vial	1 bottle	Store at 2-8°C ***for six months
Plate Cover Seals	4 pieces	

**Provided this is within the expiration date of the kit.

OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED

1. Microplate reader capable of measuring absorbance at 450 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squir bottle, manifold dispenser, or automated microplate washer.
5. 500 mL graduated cylinder.
6. Human CCL5 controls (optional; available from Solarbio).

SPECIMEN COLLECTION & STORAGE

Cell Culture Supernates - Centrifuge cell culture media at 1000×g to remove debris. Assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

Serum - Use a serum separator tube (SST) and allow samples to clot for 2 hours at room temperature or overnight at 2-8°C. Centrifuge at approximately 15 minutes at 1000×g. Assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

Plasma - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000×g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20°C. Avoid repeated freeze-thaw cycles.

REAGENTS PREPARATION

1. **Temperature returning** - Bring all kit components and specimen to room temperature (20-25°C) before use.
2. **Wash Buffer** - Dilute 30mL of Wash Buffer Concentrate with 570mL of deionized or distilled water to prepare 600mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.

3. **Standard/Specimen** - Reconstitute the Standard with 1.0mL of Standard/Sample Diluent. This reconstitution produces a stock solution of 2000 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 750uL of Standard/Sample Diluent into the 500pg/mL tube, and add 250uL stock solution of 2000 pg/mL into it to get the high standard of 500 pg/mL. Pipette 500uL of Standard/Sample Diluent into the remaining tubes. Use the high standard to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 500 pg/mL standard serves as the high standard. The Standard/Sample Diluent serves as the zero standard (0 pg/mL).


***If you do not run out of re-melting standard, store it at -20°C. Diluted standard shall not be reused.**

4. **Working solution of Biotin-Conjugate anti-human CCL5 antibody:** Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.

***The working solution should be used within one day after dilution.**

5. **Working solution of Streptavidin-HRP:** Make a 1:40 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.

***The working solution should be used within one day after dilution.**

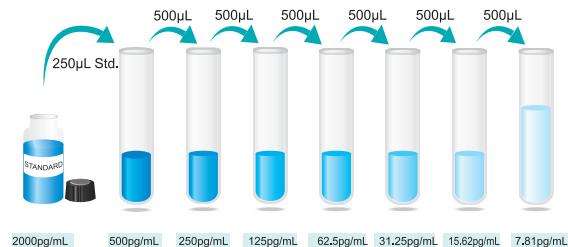

Preparation of CCL5 standard dilutions