

REFERENCES

1. Lee, J. et al.(2001). J. Biol. Chem.276:1660-1664.
2. Pan, G., D. et al. (2001). J. Immunol.167:6559-6567.
3. Hurst, S. D. et al.(2002). J. Immunol.169:443-453.
4. Aggarwal, S. et al.(2002). J. Leukocyte Biol.71:1-8.
5. Moseley, T. A. et al.(2003). Cytokine Growth Factor Rev.14:155-174.
6. Kolls, J. K. et al.(2004). Immunity 21:467-476.
7. Kim, M. R. et al.(2002). Blood 100:2330-2340.
8. Schwarzenberger, P.et al.(1998). J. Immunol. 161:6383-6389.
9. Linden, A. et al.(2000). Eur. Respir. J.15:973-977.

Mouse IL-25/IL-17EImmunoassay

Catalog Number: SEKM-0026

For the quantitative determination of mouse interleukin-25/17E(IL-25/IL-17E) concentrations in cell culture supernates, serum, and plasma.

For research use only. Not for use in diagnostic procedures.

MANUFACTURED AND DISTRIBUTED BY:

Country | Company: China | Beijing Solarbio Science & Technology Co.,Ltd
Address:NO.85A, Liandong U Valley, Tongzhou District, Beijing, P.R.China.
Tel: 86-10-56371241 Fax: 86-10-56371282 E-mail: service@solarbio.com

TABLE OF CONTENTS

SECTION	PAGE
BACKGROUND.....	01
PRINCIPLE OF THE ASSAY.....	01
TECHNICAL HINTS AND LIMITATIONS.....	02
PRECAUTIONS.....	02
KIT COMPONENTS& STORAGE CONDITIONS.....	03
OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED.....	04
SPECIMEN COLLECTION & STORAGE.....	04
REAGENTS PREPARATION.....	04
ASSAY PROCEDURE.....	06
CALCULATION OF RESULTS.....	06
PERFORMANCE CHARACTERISTICS.....	08
REFERENCES.....	10

LINEARITY: To assess the linearity of the assay, three samples were spiked with high concentrations of IL-25/IL-17E in various matrices and diluted with the appropriate Sample Diluent to produce samples with values within the dynamic range of the assay. (The plasma samples were initially diluted 1:1)

Dilution ratio	Recovery(%)	Citrate plasma	Cell culture supernatants
1:2	Average% of Expected	95	103
	Range (%)	87-102	93-106
1:4	Average% of Expected	92	94
	Range (%)	90-113	91-106
1:8	Average% of Expected	96	102
	Range (%)	89– 102	90–113
1:16	Average% of Expected	99	105
	Range (%)	90–107	94–116

Performance Characteristics

SENSITIVITY: The minimum detectable dose was 15pg/mL.

SPECIFICITY: This assay recognizes both natural and recombinant mouse IL-25/IL-17E. The factors listed below were prepared at 100ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

Factors assayed for cross-reactivity

Recombinant mouse	Recombinant rat	Recombinant human
IL-17		
IL-17B		
IL-17C		
IL-17D		
IL-17F		

REPEATABILITY: The coefficient of variation of both intra-assay and inter-assay were less than 10%.

RECOVERY: The recovery of IL-25/IL-17E spiked to three different levels in four samples throughout the range of the assay in various matrices was evaluated.

Recovery of IL-25/IL-17E in two matrices

Sample Type	Average % of Expected Range (%)	Range (%)
Citrate plasma	95	91-110
Cell culture supernatants	103	95-103

BACKGROUND


The IL-17 family is comprised of at least six proinflammatory cytokines that share a conserved cysteine-knot structure but diverge at the N-terminus. IL-17 family members are glycoproteins secreted as dimers that induce local cytokine production and recruit granulocytes to sites of inflammation. IL-17 is induced by IL-15 and IL-23, mainly in activated CD4+ T cells distinct from Th1 or Th2 cells. IL-17F is the most homologous to IL-17, but is induced only by IL-23 in activated monocytes. IL-17B binds the IL-17B receptor, but not the IL-17 receptor; it is most homologous with IL-17D, which is expressed by resting CD4+ T cells and CD19+ B cells. IL-17E is mainly produced by Th2 cells and recruits eosinophils to lung tissue. IL-17C has a very restricted expression pattern but has been detected in adult prostate and fetal kidney libraries.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for IL-25/IL-17E has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any IL-25/IL-17E present is captured by the coated antibody after incubation. Following extensive washing, a biotin-conjugate antibody specific for IL-25/IL-17E is added to detect the captured IL-25/IL-17E protein in sample. For signal development, horseradish peroxidase (HRP)-conjugated Streptavidin is added, followed by tetramethyl-benzidine (TMB) reagent. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Solution containing sulfuric acid is used to stop color development and the color intensity which is proportional to the quantity of bound protein is measurable at 450nm.

DESCRIPTION

Schematic diagram:

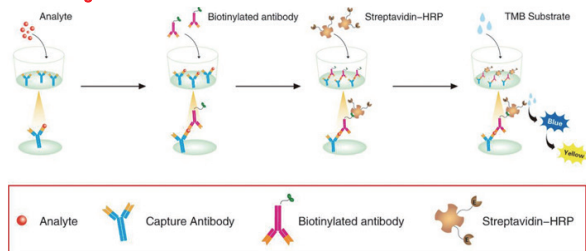

TECHNICAL HINTS AND LIMITATIONS

1. This Solarbio ELISA should not be used beyond the expiration data on the kit label.
2. To avoid cross-contamination, use a fresh reagent reservoir and pipette tips for each step.
3. To ensure accurate results, some details, such as technique, plasticware and water sources should be emphasized.
4. A thorough and consistent wash technique is essential for proper assay performance.
5. A standard curve should be generated for each set of samples assayed.
6. It is recommended that all standards and samples be assayed in duplicate.
7. Avoid microbial contamination of reagents and buffers. Buffers containing protein should be made under aseptic conditions and be prepared fresh daily.
8. In order to ensure the accuracy of the results, the standard curve should be made every time.

PRECAUTIONS


The Stop Solution suggested for use with this kit is an acid solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

DESCRIPTION

- precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
5. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Typical data using the IL-25/IL-17E ELISA

Standard(pg/ml)	OD.	OD.	Average	Corrected
0	0.053	0.059	0.056	-
31.25	0.118	0.109	0.113	0.057
62.5	0.215	0.235	0.225	0.169
125	0.361	0.389	0.375	0.319
250	0.578	0.542	0.560	0.504
500	0.896	0.902	0.899	0.843
1000	1.521	1.546	1.533	1.477
2000	2.631	2.405	2.518	2.462

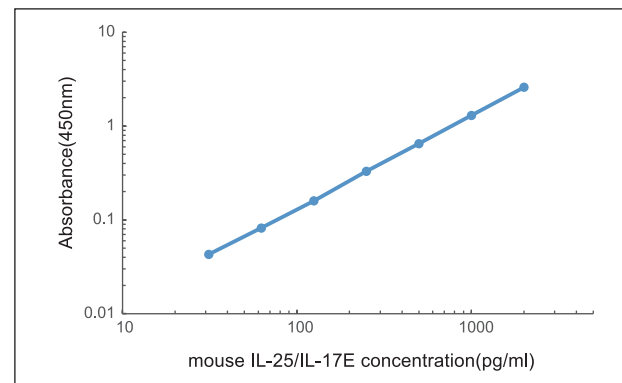

Representative standard curve for IL-25/IL-17E ELISA.

ASSAY PROCEDURE

Prepare all reagents and standards as directed. Wash the plate 3 times before assay.



Add 100µl standard or samples to each well, incubate 90 minutes, 37°C.



Aspirate and wash 4 times

Add 100µl working solution of Biotin-Conjugate anti-mouse IL-25/IL-17E antibody to each well, incubate 60 minutes, 37°C.



Aspirate and wash 4 times

Add 100µl working solution of Streptavidin-HRP to each well, incubate 30 minutes, 37°C.



Aspirate and wash 5 times

Add 100µl Substrate solution to each well, incubate 15 minutes, 37°C.
Protect from light.



Add 50µl Stop solution to each well. Read at 450nm within 5 minutes.

CALCULATION OF RESULTS

1. The standard curve is used to determine the amount of specimens.
2. First, average the duplicate readings for each standard, control, and sample. All O.D. values are subtracted by the mean value of blank control before result interpretation.
3. Construct a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
4. The data may be linearized by plotting the log of the IL-25/IL-17E concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less

KIT COMPONENTS & STORAGE CONDITIONS

PART	SIZE	STORAGE OF OPENED/ RECONSTITUTED MATERIAL
Microwell Plate-antibody coated 96-well Microplate (8 wells x12 strips)	1 plate	Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2 – 8°C**
Standard - lyophilized, 2000 pg/ml upon reconstitution	2 vials	Aliquot and Store at -20°C** for six months
Concentrated Biotin-Conjugated antibody (100X) - 120 ul/vial	1 vial	Store at 2-8°C **for six months
Concentrated Streptavidin-HRP solution (100X) - 120 ul/vial	1 vial	Store at 2-8°C **for six months
Standard /sample Diluent - 16ml/vial	1 bottle	Store at 2-8°C **for six months
Biotin-Conjugate antibody Diluent- 16ml/vial	1 bottle	Store at 2-8°C **for six months
Streptavidin-HRP Diluent - 16ml/vial	1 bottle	Store at 2-8°C **for six months
Wash Buffer Concentrate (20x) - 30 ml/vial	1 bottle	Store at 2-8°C **for six months
Substrate Solution - 12 ml/vial	1 bottle	Store at 2-8°C **for six months
Stop Solution - 12 ml/vial	1 bottle	Store at 2-8°C **for six months
Plate Cover Seals	4 pieces	

**Provided this is within the expiration date of the kit.

OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED

1. Microplate reader capable of measuring absorbance at 450 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squirrt bottle, manifold dispenser, or automated microplate washer.
5. 500 mL graduated cylinder.

SPECIMEN COLLECTION & STORAGE

Cell Culture Supernates- Centrifuge cell culture media at 1000×g to remove debris. Assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.

Serum - Use a serum separator tube (SST) and allow samples to clot for 2 hours at room temperature or overnight at 2–8°C. Centrifuge approximately for 15 minutes at 1000×g. Assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.


Plasma - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000×g within 30 minutes of collection. Assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.

Note: It is recommended to conduct a pre-test before the formal experiment to determine the dilution ratio.

REAGENTS PREPARATION

1. **Temperature returning** - Bring all kit components and specimen to room temperature (20–25°C) before use.
2. **Wash Buffer** - Dilute 30mL of Wash Buffer Concentrate with 570mL of deionized or distilled water to prepare 600mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.

3. **Standard/Sample**- Reconstitute the Standard with 1.0mL of Standard/Sample Diluent. This reconstitution produces a stock solution of 2000pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 500µL of Standard/Sample Diluent into 1000pg/mL tube and the remaining tubes. Use the stock solution of 2000pg/mL to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 2000pg/mL standard serves as the high standard. The Standard/specimen Diluent serves as the zero standard (0 pg/mL).

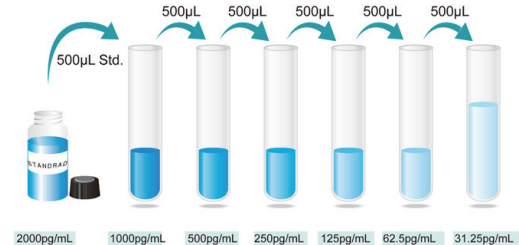**Preparation of IL-25/IL-17E standard dilutions**

***If you do not run out of re-melting standard, store it at -20°C . Diluted standard shall not be reused.**

4. **Working solution of Biotin-Conjugate anti-mouse IL-25/IL-17E antibody:** Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.
***The working solution should be used within one day after dilution.**
5. **Working solution of Streptavidin-HRP:** Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.
***The working solution should be used within one day after dilution.**