

**REFERENCES**

- 1.Ealick, S.E. et al. (1991) Science 252:698.
- 2.Rashidbaigi, A. et al. (1986) Proc. Natl. Acad. Sci. USA 83:384.
- 3.Pfizenmaier, K. et al. (1988) J. Immunol. 141:856.
- 4.Schreiber, R.D. and M. Aguet (1994) in Guidebook to Cytokines and their Receptors, N.A. Nicola ed., Oxford University Press, New York, p. 120.
- 5.Aguet, M. et al. (1988) Cell 55:273.
- 6.Weil, J. et al. (1983) Nature 301:437.
- 7.Harris, C.A. et al. (1992) J. Biol. Chem. 267:6865.
- 8.Williams, B.R.G. (1991) Eur. J. Biochem. 200:1.
- 9.Shuai, K. et al. (1993) Science 261:1744.

## Mouse IFN- $\gamma$ Immunoassay

Catalog Number: SEKM-0031

For the quantitative determination of mouse Interferon-gamma (IFN- $\gamma$ ) concentrations in cell culture supernates, serum, and plasma.

For research use only. Not for use in diagnostic procedures.

### MANUFACTURED AND DISTRIBUTED BY:

Country | Company: China | Beijing Solarbio Science & Technology Co.,Ltd  
Address:NO.85A, Liandong U Valley, Tongzhou District, Beijing, P.R.China.  
Tel: 86-010-50973105 Fax: 86-10-56371282 E-mail: service@solarbio.com

## TABLE OF CONTENTS

SECTION	PAGE
BACKGROUND.....	01
PRINCIPLE OF THE ASSAY.....	01
TECHNICAL HINTS AND LIMITATIONS.....	02
PRECAUTIONS.....	02
KIT COMPONENTS& STORAGE CONDITIONS.....	03
OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED.....	04
SPECIMEN COLLECTION & STORAGE.....	04
REAGENTS PREPARATION.....	04
ASSAY PROCEDURE.....	06
CALCULATION OF RESULTS.....	06
PERFORMANCE CHARACTERISTICS.....	08
REFERENCES.....	10

**LINEARITY:** To assess the linearity of the assay, three samples were spiked with high concentrations of IFN- $\gamma$  in various matrices and diluted with the appropriate Sample Diluent to produce samples with values within the dynamic range of the assay.

Dilution ratio	Recovery(%)	Citrate plasma	Cell culture supernatants
1:2	Average% of Expected	93	95
	Range (%)	85-98	90-112
1:4	Average% of Expected	92	97
	Range (%)	86-103	91-105
1:8	Average% of Expected	105	99
	Range (%)	95–114	95–103
1:16	Average% of Expected	107	104
	Range (%)	98–116	100–107

**Performance Characteristics**

**SENSITIVITY:** The minimum detectable dose was 3.9 pg/mL.

**SPECIFICITY:** This assay recognizes both natural and recombinant mouse IFN- $\gamma$ . The factors listed below were prepared at 100ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

Factors assayed for cross-reactivity

Recombinant mouse	Recombinant rat	Recombinant human
IL-2		IFN- $\gamma$
IL-4		
IL-5		
IL-6		
IL-7		
IL-9		
IL-1 $\alpha$		
SCF		
VEGF		

**REPEATABILITY:** The coefficient of variation of both intra-assay and inter-assay were less than 10%.

**RECOVERY:** The recovery of IFN- $\gamma$  spiked to three different levels in four samples throughout the range of the assay in various matrices was evaluated.

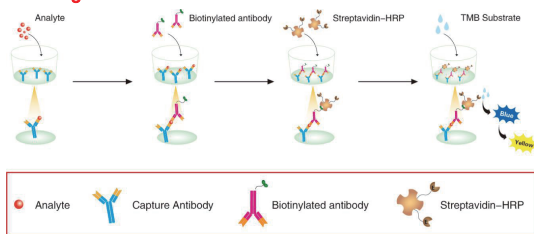
Sample Type	Average % of Expected Range(%)	Range(%)
Citrate plasma	89	82-103
Cell culture supernatants	93	90-113

**BACKGROUND**

IFN-gamma (Interferon-gamma) is the prototype proinflammatory cytokine and is produced by a variety of immune cells under inflammatory conditions, notably by T cells and NK cells. It plays a key role in host defense by promoting the development and activation of Th1 cells, chemoattraction and activation of monocytes and macrophages, upregulation of antigen presentation molecules, and immunoglobulin class switching in B cells. It also exhibits antiviral, antiproliferative, and apoptotic effects. In addition, IFN-gamma functions as an anti-inflammatory mediator by promoting the development of regulatory T cells and inhibiting Th17 cell differentiation. IFN-gamma dimers signal through a receptor complex of two IFN-gamma R1 and two IFN-gamma R2 subunits.

**PRINCIPLE OF THE ASSAY**

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for IFN- $\gamma$  has been pre-coated onto a microplate. Standards and samples are pipetted into the wells; any IFN- $\gamma$  present is captured by the coated antibody after incubation. After washing away any unbound substances, a biotin-conjugate antibody specific for IFN- $\gamma$  is added to detect the captured IFN- $\gamma$  protein in the sample. Following a wash to remove any unbound combination, horseradish peroxidase (HRP)-conjugated Streptavidin is added to the wells. After extensive washing, a tetramethyl-benzidine (TMB) reagent is added to the wells for signal development. Solution containing sulfuric acid is used to stop color development. The color intensity, proportional to the quantity of bound protein, is then measurable at 450nm.

**Schematic diagram:****TECHNICAL HINTS AND LIMITATIONS**

1. This Solarbio ELISA should not be used beyond the expiration data on the kit label.
2. To avoid cross-contamination, use a fresh reagent reservoir and pipette tips for each step.
3. To ensure accurate results, some details, such as technique, plasticware and water sources should be emphasized.
4. A thorough and consistent wash technique is essential for proper assay performance.
5. A standard curve should be generated for each set of samples assayed.
6. It is recommended that all standards and samples be assayed in duplicate.
7. Avoid microbial contamination of reagents and buffers. Buffers containing protein should be made under aseptic conditions and be prepared fresh daily.
8. In order to ensure the accuracy of the results, the standard curve should be made every time.

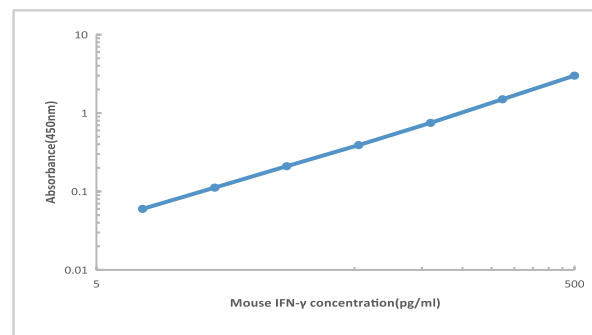
**PRECAUTIONS**

The Stop Solution suggested for use with this kit is an acid solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

4. The data may be linearized by plotting the log of the IFN- $\gamma$  concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
5. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Typical data using the IFN- $\gamma$  ELISA

Standardized (pg/ml)	OD.	OD.	Average	Corrected
0	0.105	0.100	0.103	-
7.81	0.161	0.168	0.165	0.062
15.63	0.235	0.244	0.240	0.137
31.25	0.349	0.363	0.356	0.254
62.5	0.554	0.577	0.566	0.463
125	0.901	0.937	0.919	0.817
250	1.466	1.526	1.496	1.394
500	2.389	2.487	2.438	2.336

Representative standard curve for IFN- $\gamma$  ELISA.

**ASSAY PROCEDURE**

Prepare all reagents and standards as directed. Wash the plate 3 times before assay.



Add 100µl standard or samples to each well, incubate 90 minutes, 37°C.



Aspirate and wash 4 times

Add 100µl working solution of Biotin-Conjugate anti-mouse IFN-γ antibody to each well, incubate 60 minutes, 37°C.



Aspirate and wash 4 times

Add 100µl working solution of Streptavidin-HRP to each well, incubate 30 minutes, 37°C.



Aspirate and wash 5 times

Add 100µl Substrate solution to each well, incubate 15 minutes, 37°C. Protect from light.



Add 50µl Stop solution to each well. Read at 450nm within 5 minutes.

**CALCULATION OF RESULTS**

1. The standard curve is used to determine the amount of specimens.
2. First, average the duplicate readings for each standard, control, and sample. All O.D. values are subtracted by the mean value of blank control before result interpretation.
3. Construct a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.

**KIT COMPONENTS & STORAGE CONDITIONS**

PART	SIZE	STORAGE OF OPENED/ RECONSTITUTED MATERIAL
Microwell Plate - antibody coated 96-well Microplate (8 wells x12 strips)	1 plate	Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2 – 8°C**
Standard - lyophilized, 2000 pg/ml upon reconstitution	2 vials	Store at 2-8°C ***for six months
Concentrated Biotin-Conjugated antibody(100X) - 120 µl/vial	1 vial	Store at 2-8°C ***for six months
Concentrated Streptavidin-HRP solution(100X) - 120 µl/vial	1 vial	Store at 2-8°C ***for six months
Standard /sample Diluent - 16 ml/vial	1 bottle	Store at 2-8°C ***for six months
Biotin-Conjugate antibody Diluent - 16 ml/vial	1 bottle	Store at 2-8°C ***for six months
Streptavidin-HRP Diluent - 16 ml/vial	1 bottle	Store at 2-8°C ***for six months
Wash Buffer Concentrate (20x) - 30 ml/vial	1 bottle	Store at 2-8°C ***for six months
Substrate Solution - 12 ml/vial	1 bottle	Store at 2-8°C ***for six months
Stop Solution - 12 ml/vial	1 bottle	Store at 2-8°C ***for six months
Plate Cover Seals	4 bottle	

\*\*Provided this is within the expiration date of the kit.

**OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED**

1. Microplate reader capable of measuring absorbance at 450 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squirrt bottle, manifold dispenser, or automated microplate washer.
5. 500 mL graduated cylinder.

**SPECIMEN COLLECTION & STORAGE**

**Cell Culture Supernates** - Centrifuge cell culture media at 1000×g to remove debris. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

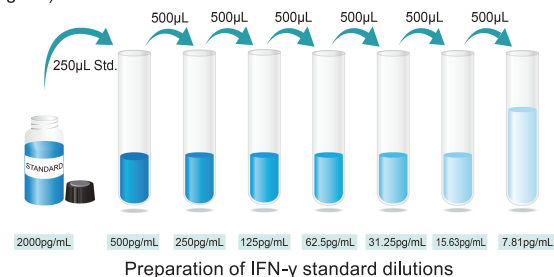
**Serum** - Use a serum separator tube (SST) and allow samples to clot for 2 hours at room temperature or overnight at 2-8°C. Centrifuge approximately for 15 minutes at 1000×g. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

**Plasma** - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000×g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

**REAGENTS PREPARATION**

1. **Temperature returning** - Bring all kit components and specimen to room temperature (20-25°C) before use.
2. **Wash Buffer** - Dilute 30mL of Wash Buffer Concentrate with 570mL of deionized or distilled water to prepare 600mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.

3. **Standard/Sample** - Reconstitute the Standard with 1.0mL of Standard/Sample Diluent. This reconstitution produces a stock solution of 2000pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 750µL of Standard/Sample Diluent into the 500pg/mL tube, and add 250µL stock solution of 2000pg/mL into it to get the 500pg/mL. Use the 500pg/mL to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 500 pg/mL standard serves as the high standard. The Standard/Sample Diluent serves as the zero standard (0 pg/mL).



**\*If you do not run out of re-melting standard, store it at -20°C. Diluted standard shall not be reused.**

4. **Working solution of Biotin-Conjugate anti-mouse IFN-γ antibody:** Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.
5. **Working solution of Streptavidin-HRP:** Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.