

LINEARITY: To assess the linearity of the assay, three samples were spiked with high concentrations of CCL28 in various matrices and diluted with the appropriate Sample Diluent to produce samples with values within the dynamic range of the assay. (The plasma samples were initially diluted 1:1)

Dilution ratio	Recovery (%)	Citrate plasma	Cell culture supernatants
1:2	Average% of Expected	96	98
	Range (%)	90–101	92–104
1:4	Average% of Expected	98	101
	Range (%)	91–105	94–108
1:8	Average% of Expected	101	96
	Range (%)	95–107	88–104
1:16	Average% of Expected	104	99
	Range (%)	99–109	92–106

REFERENCES

1.Wang, W. et al. (2000) J. Biol. Chem. 275:22313.


Mouse CCL28 Immunoassay

Catalog Number: SEKM-0045

For the quantitative determination of mouse CCL28 concentrations in cell culture supernates, serum, and plasma.

For research use only. Not for use in diagnostic procedures.

Country | Company: China | Beijing Solarbio Science & Technology Co.,Ltd
 Address:NO.85A, Liandong U Valley, Tongzhou District, Beijing, P.R.China.
 Tel: 86-10-56371241 Fax: 86-10-56371282 E-mail: service@solarbio.com

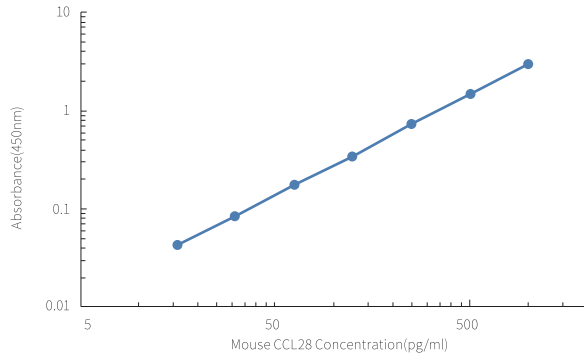

Representative standard curve for CCL28 ELISA.

BACKGROUND

Mouse CCL28 (CC chemokine ligand 28) is a novel CC chemokine cloned from a Rag-1 mouse kidney cDNA library. Mouse CCL28 cDNA encodes a 130 amino acid (aa) residue precursor protein with a putative 22 aa residue signal peptide that is cleaved to produce the 108 aa residue mature protein. Human and mouse CCL28 are highly conserved, sharing 83% aa identity in their mature regions. Among CC chemokines, CCL28 shares the most homology with CCL27/CTACK. The mouse CCL28 gene has been mapped to the distal region of chromosome 13. Mouse CCL28 is produced by epithelial cells. Based on Northern blot analysis, it is mainly expressed in testes, kidney and brain. The receptor for CCL28 has been identified as the CCR10 which is also the receptor for CCL27/CTACK.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for CCL28 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any CCL28 present is captured by the coated antibody after incubation. Following extensive washing, a biotin-conjugate antibody specific for CCL28 is added to detect the captured CCL28 protein in sample. For signal development, horseradish peroxidase (HRP)-conjugated Streptavidin is added, followed by tetramethyl-benzidine (TMB) reagent. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Solution containing sulfuric acid is used to stop color development and the color intensity which is proportional to the quantity of bound protein is measurable at 450nm.



TECHNICAL HINTS AND LIMITATIONS

1. This Solarbio ELISA should not be used beyond the expiration date on the kit label.
2. To avoid cross-contamination, use a fresh reagent reservoir and pipette tips for each step.
3. To ensure accurate results, some details, such as technique, plasticware and water sources should be emphasized.
4. A thorough and consistent wash technique is essential for proper assay performance.
5. A standard curve should be generated for each set of samples assayed.
6. It is recommended that all standards and samples be assayed in duplicate.
7. Avoid microbial contamination of reagents and buffers. Buffers containing protein should be made under aseptic conditions and be prepared fresh daily.
8. In order to ensure the accuracy of the results, the standard curve should be made every time.

PRECAUTIONS

The Stop Solution suggested for use with this kit is an acid solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

CALCULATION OF RESULTS

1. The standard curve is used to determine the amount of specimens.
2. First, average the duplicate readings for each standard, control, and sample. All O.D. values are subtracted by the mean value of blank control before result interpretation.
3. Construct a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
4. The data may be linearized by plotting the log of the CCL28 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
5. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Standard (pg/ml)	OD.	OD.	Average	Corrected
0	0.018	0.022	0.020	—
15.625	0.250	0.253	0.252	0.232
31.25	0.289	0.292	0.290	0.270
62.5	0.429	0.434	0.432	0.412
125	0.682	0.690	0.686	0.666
250	1.108	1.120	1.114	1.094
500	1.803	1.823	1.813	1.793
1000	2.938	2.971	2.955	2.935

5. Working solution of Biotin-Conjugate anti-mouse CCL28 antibody:

Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.

***The working solution should be used within one day after dilution.**

ASSAY PROCEDURE

Prepare all reagents and standards as directed. Wash the plate 3 times before assay.



Add 100µl standard or samples to each well, incubate 90 minutes, 37°C.



Aspirate and wash 4 times

Add 100µl working solution of Biotin-Conjugate anti-mouse CCL28 antibody to each well, incubate 60 minutes, 37°C.



Aspirate and wash 4 times

Add 100µl working solution of Streptavidin-HRP to each well, incubate 30 minutes, 37°C.



Aspirate and wash 5 times

Add 100µl Substrate solution to each well, incubate 15 minutes, 37°C.
Protect from light.



Add 50µl Stop solution to each well. Read at 450nm within 5 minutes.

KIT COMPONENTS & STORAGE CONDITIONS

PART	SIZE	STORAGE OF OPENED/ RECONSTITUTED MATERIAL
Microwell Plate - antibody coated 96-well Microplate (8 wells ×12 strips)	1 plate	Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2 – 8°C**
Standard -lyophilized, 2000 pg/vial upon reconstitution	2 vials	Aliquot and Store at -20°C** for six months
Concentrated Biotin-Conjugated antibody (100X) - 120 ul/vial	1 vial	Store at 2-8°C ** for six months
Concentrated Streptavidin-HRP solution (100X) - 120 ul/vial	1 vial	Store at 2-8°C** for six months
Standard /Sample Diluent - 16 ml/vial	1 bottle	Store at 2-8°C** for six months
Biotin-Conjugate antibody Diluent - 16 ml/vial	1 bottle	Store at 2-8°C** for six months
Streptavidin-HRP Diluent - 16 ml/vial	1 bottle	Store at 2-8°C** for six months
Wash Buffer Concentrate (20x) - 30 ml/vial	1 bottle	Store at 2-8°C** for six months
Substrate Solution - 12 ml/vial	1 bottle	Store at 2-8°C** for six months
Stop Solution - 12 ml/vial	1 bottle	Store at 2-8°C** for six months
Plate Cover Seals	4 pieces	

**Provided this is within the expiration date of the kit.

OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED

1. Microplate reader capable of measuring absorbance at 450 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squir bottle, manifold dispenser, or automated microplate washer.
5. 500 mL graduated cylinder.

SPECIMEN COLLECTION & STORAGE

Cell Culture Supernates - Centrifuge cell culture media at 1000×g to remove debris. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Serum - Use a serum separator tube (SST) and allow samples to clot for 2 hours at room temperature or overnight at 2-8°C. Centrifuge approximately for 15 minutes at 1000×g. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Plasma - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000×g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Note: The normal mouse serum or plasma samples are suggested to make a 1:2 dilution.

REAGENTS PREPARATION


5. **Working solution of Streptavidin-HRP: Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.**

***The working solution should be used within one day after dilution.**

Temperature returning - Bring all kit components and specimen to room temperature (20-25°C) before use.

3. **Wash Buffer** - Dilute 30mL of Wash Buffer Concentrate with 570mL of deionized or distilled water to prepare 600mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.

Standard/Specimen - Reconstitute the Standard with 1.0mL of Standard/Sample Diluent. This reconstitution produces a stock solution of 2000 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 500µL of Standard/Sample Diluent into 1000pg/ml tube and the remaining tubes. Use the

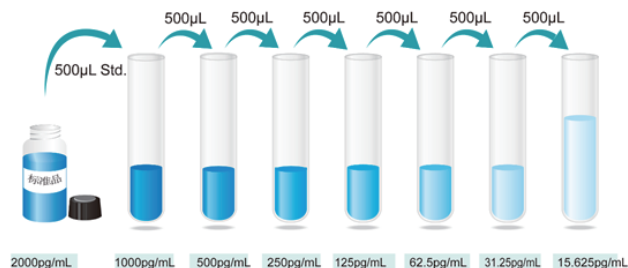

- stock solution of 1000pg/mL to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 1000 pg/mL standard serves as the high standard. The Standard/
4. Sample Diluent serves as the zero standard (0 pg/mL).

Preparation of CCL28 standard dilutions

*If you do not run out of re-melting standard, store it at -20°C. Diluted standard shall not be reused.