

**REFERENCES**

- 1.Park, A.Y. and P. Scott (2001) Scand. J. Immunol. 53:529.
2. Trinchieri, G. et al. (2003) Immunity 19:641.
3. Brombacher, F. et al. (2003) Trends Immunol. 24:207.
4. Lankford, C.S. and D.M. Frucht (2003) J. Leukoc. Biol. 73:49
5. Temblay JN, et al. J Allergy Clin Immunol. 2007 Sep;120(3):

## Rabbit IL-12 Immunoassay

Catalog Number: SEKRT-0404

For the quantitative determination of Rabbit IL-12 concentrations in cell culture supernates, serum, and plasma.

For research use only. Not for use in diagnostic procedures.

Country | Company: China | Beijing Solarbio Science & Technology Co.,Ltd  
Address:NO.85A, Liandong U Valley, Tongzhou District, Beijing, P.R.China.  
Tel: 86-10-56371241 Fax: 86-10-56371282 E-mail: [service@solarbio.com](mailto:service@solarbio.com)

## TABLE OF CONTENTS

| SECTION                                       | PAGE |
|---|------|
| BACKGROUND.....                               | 01   |
| PRINCIPLE OF THE ASSAY.....                   | 01   |
| TECHNICAL HINTS AND LIMITATIONS.....          | 02   |
| PRECAUTIONS.....                              | 02   |
| KIT COMPONENTS& STORAGE CONDITIONS.....       | 03   |
| OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED..... | 04   |
| SPECIMEN COLLECTION & STORAGE.....            | 04   |
| REAGENTS PREPARATION.....                     | 04   |
| ASSAY PROCEDURE.....                          | 06   |
| CALCULATION OF RESULTS.....                   | 07   |
| PERFORMANCE CHARACTERISTICS.....              | 08   |
| REFERENCES.....                               | 10   |

**REPEATABILITY:** The coefficient of variation of both intra-assay and inter-assay were less than 10%.


**RECOVERY:** The recovery of IL-12 spiked to three different levels in four samples throughout the range of the assay in various matrices was evaluated.

### Recovery of IL-12 in two matrices

| Sample Type               | Average % of Expected Range (%) | Range (%) |
|---------------------------|---------------------------------|-----------|
| Citrate plasma            | 95                              | 88-105    |
| Cell culture supernatants | 97                              | 86-103    |

**LINEARITY:** To assess the linearity of the assay, three samples were spiked with high concentrations of IL-12 in various matrices and diluted with the appropriate Sample Diluent to produce samples with values within the dynamic range of the assay. (The plasma samples were initially diluted 1:1)

| Dilution ratio | Recovery(%)          | Citrate plasma | Cell culture supernatants |
|----------------|----------------------|----------------|---------------------------|
| 1:2            | Average% of Expected | 92             | 107                       |
|                | Range (%)            | 85-101         | 99-114                    |
| 1:4            | Average% of Expected | 93             | 109                       |
|                | Range (%)            | 83-106         | 101-120                   |

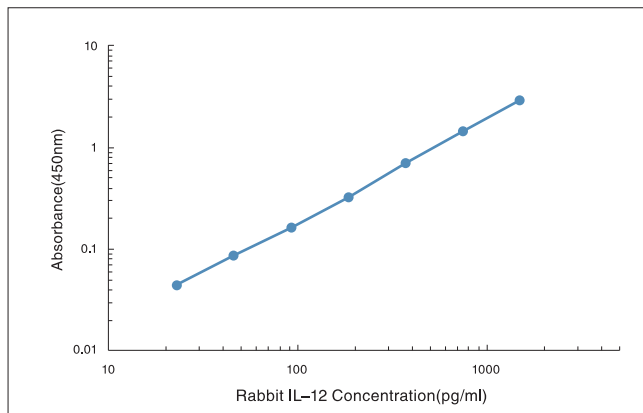

Representative standard curve for IL-12 ELISA

### Performance Characteristics

**SENSITIVITY:**The minimum detectable dose was 5 pg/mL.

**SPECIFICITY:**This assay recognizes both natural and recombinant Rabbit IL-12. The factors listed below were prepared at 10ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

BMP1,BMP2,BMP4,HGF,IL-1 $\beta$ ,IL-2,IL-4,IL-5,IL-6,IL-8,IL-10,IL-13,IL-15, IFN,MMP-2,TGF $\beta$ 1,TGF $\beta$ 2,TGF $\beta$ 3,TLR1,TLR2,TLR3,TNF- $\alpha$ ,VEGF

### BACKGROUND

Interleukin 12 (IL12) is the founding member of the IL12 family of heterodimeric cytokines, which have important immunological functions. IL12 is composed of two disulfidelinked subunits of 35 kDa and 40 kDa, respectively. The 35 kDa subunit (p35) is an helical protein homologous to IL6 and GCSF. The 40 kDa subunit (p40) contains one fibronectin type III and one Ig C2like domain and has a high degree of structural homology to type I cytokine receptors. Whereas p35 subunit is unique to IL 12, the p40 subunit is also utilized in IL23. Mature rat p40 contains 313 aa and can exist in multiple forms, including monomer, homodimer, heterodimer linked to p19 (forming IL 23), and heterodimer linked to p35 (forming IL12). The expression of p40 is upregulated by substances such as LPS and CpG that activate antigen-presenting cells.


Cells known to produce IL12 include macrophages, dendritic cells, monocytes, Langerhans cells, neutrophils, and keratinocytes. The activities of IL12 are mediated by the receptor complex composed of two type I transmembrane proteins: a 100 kDa ligandbinding subunit (IL12 R $\beta$ 1) and a 130 kDa signal transducing subunit (IL12 R $\beta$ 2). IL-12 mediates enhancement of the cytotoxic activity of NK cells and CD8+ cytotoxic T lymphocytes and potentiates the expansion and late activation of Th1 CD4+ T cells. IL-12 also has anti-angiogenic activity . IL-12 is linked with autoimmunity, and normally stops allergies to food developing.

### PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for IL-12 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any IL-12 present is captured by the coated antibody after incubation. Following extensive washing, a biotin-conjugate antibody specific for IL-12 is added to detect the captured IL-12 protein in sample. For signal development, horseradish peroxidase (HRP)-conjugated Streptavidin is added, followed by tetramethyl-benzidine (TMB) reagent. Following a wash to remove any

## DESCRIPTION

unbound combination, and enzyme conjugate is added to the wells. Solution containing sulfuric acid is used to stop color development and the color intensity which is proportional to the quantity of bound protein is measurable at 450nm.

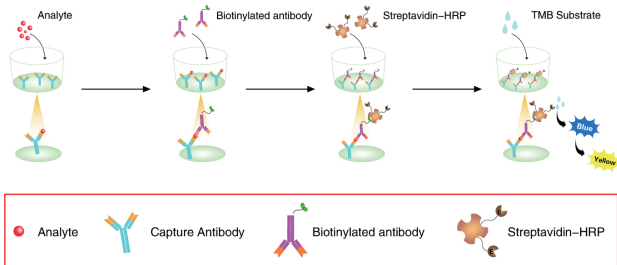

## TECHNICAL HINTS AND LIMITATIONS

- 1.This Solarbio ELISA should not be used beyond the expiration data on the kit label.
- 2.To avoid cross-contamination, use a fresh reagent reservoir and pipette tips for each step.
- 3.To ensure accurate results, some details, such as technique, plasticware and water sources should be emphasized.
- 4.A thorough and consistent wash technique is essential for proper assay performance.
- 5.A standard curve should be generated for each set of samples assayed.
- 6.It is recommended that all standards and samples be assayed in duplicate.
- 7.Avoid microbial contamination of reagents and buffers. Buffers containing protein should be made under aseptic conditions and be prepared fresh daily.
- 8.In order to ensure the accuracy of the results, the standard curve should be made every time.

## DESCRIPTION

### CALCULATION OF RESULTS

1. The standard curve is used to determine the amount of specimens.
2. First, average the duplicate readings for each standard, control, and sample. All O.D. values are subtracted by the mean value of blank control before result interpretation.
3. Construct a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
4. The data may be linearized by plotting the log of the IL-12 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
5. This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Typical data using the IL-12 ELISA

| Std (pg/mL) | O.D.1 | O.D.2 | Averag | Correct |
|-------------|-------|-------|--------|---------|
| 0           | 0.061 | 0.063 | 0.062  | ---     |
| 25          | 0.092 | 0.094 | 0.093  | 0.031   |
| 50          | 0.175 | 0.183 | 0.179  | 0.117   |
| 100         | 0.282 | 0.275 | 0.278  | 0.216   |
| 200         | 0.459 | 0.471 | 0.465  | 0.403   |
| 400         | 0.863 | 0.882 | 0.872  | 0.810   |
| 800         | 1.487 | 1.459 | 1.473  | 1.411   |
| 1600        | 2.546 | 2.570 | 2.558  | 2.496   |

## DESCRIPTION

Antibody and mix thoroughly prior to the assay. If the partial antibody is used, make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.

**\*The working solution should be used within one day after dilution.**

5. **Working solution of Streptavidin-HRP(120µL)**–Centrifuge for 1 min at 6000 x g to bring down the material prior to open the vial. The vial contains 120µL HRP Conjugate sufficient for a 96-well plate. Make 1:100 dilutions in Reagent Diluent. If the entire 96-well plate is used, add 100 µl of HRP Conjugate to 10mL of Streptavidin-HRP Diluent to make working dilution of HRP Conjugate and mix thoroughly prior to the assay. The rest of undiluted HRP Conjugate can be stored at 4°C for up to 6 months. DO NOT FREEZE.

**\*The working solution should be used within one day after dilution.**

## ASSAY PROCEDURE

Prepare all reagents and standards as directed. Wash the plate 3 times before assay.



Add 100µl standard or samples to each well, shaking with Micro-oscillator (100r/min) to incubate 60 minutes at room temperature(25±2°C).



Aspirate and wash 4 times

Add 100µl working solution of Biotin-Conjugate anti-Rabbit IL-12 antibody to each well, shaking with Micro-oscillator (100r/min) to incubate 60 minutes at room temperature(25±2°C).



Aspirate and wash 4 times

Add 100µl working solution of Streptavidin-HRP to each well, shaking with Micro-oscillator (100r/min) to incubate 20 minutes at room temperature(25±2°C).



Aspirate and wash 5 times

Add 100µl Substrate solution to each well, incubate 5-20 minutes (depending on signal) at room temperature(25±2°C). Protect from light.



Add 50µl Stop solution to each well. Read at 450nm within 5 minutes.

## DESCRIPTION

### PRECAUTIONS

The Stop Solution suggested for use with this kit is an acid solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

### KIT COMPONENTS & STORAGE CONDITIONS

| PART  | SIZE     | STORAGE OF OPENED/ RECONSTITUTED MATERIAL   |
|---|----------|---|
| Microwell Plate - antibody coated 96-well Microplate (8 wells x12 strips) | 1 plate  | Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2 – 8 C** |
| Standard - lyophilized, 1600pg/ml upon reconstitution                     | 2 vials  | Aliquot and Store at -20°C** for six months   |
| lyophilized Biotin-Conjugated antibody                                    | 1 vial   | Store at 2-8°C **for six months   |
| Concentrated Streptavidin-HRP   | 1 vial   | Store at 2-8°C **for six months   |
| Standard /sample Diluent  | 1 bottle | Store at 2-8°C **for six months   |
| Biotin-Conjugate antibody Diluent   | 1 bottle | Store at 2-8°C **for six months   |
| Streptavidin-HRP Diluent  | 1 bottle | Store at 2-8°C **for six months   |
| 20 x Wash Buffer Concentrate  | 1 bottle | Store at 2-8°C **for six months   |
| Substrate Solution  | 1 bottle | Store at 2-8°C **for six months   |
| Stop Solution   | 1 bottle | Store at 2-8°C **for six months   |
| Plate Cover Seals   | 4 pieces |   |

\*\*Provided this is within the expiration date of the kit.

**OTHER SUPPLIES REQUIRED BUT NOT SUPPLIED**

1. Microplate reader capable of measuring absorbance at 450 nm.
2. Pipettes and pipette tips.
3. Deionized or distilled water.
4. Squirrt bottle, manifold dispenser, or automated microplate washer.
5. 500 mL graduated cylinder.

**SPECIMEN COLLECTION & STORAGE**

**Cell Culture Supernates** - Centrifuge cell culture media at 1000×g to remove debris. Assay immediately or aliquot and store samples at  $\leq -20^{\circ}\text{C}$ . Avoid repeated freeze-thaw cycles.

**Serum** - Use a serum separator tube (SST) and allow samples to clot for 2 hours at room temperature or overnight at 2-8°C. Centrifuge at approximately for 15 minutes at 1000×g. Assay immediately or aliquot and store samples at  $\leq -20^{\circ}\text{C}$ . Avoid repeated freeze-thaw cycles.


**Plasma** - Collect plasma using EDTA, heparin, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000×g within 30 minutes of collection. Assay immediately or aliquot and store samples at  $\leq -20^{\circ}\text{C}$ . Avoid repeated freeze-thaw cycles.

**Note: The normal Rabbit serum or plasma samples are suggested to make a 1:2 dilution.**

**REAGENTS PREPARATION**

1. **Temperature returning** - Bring all kit components and specimen to room temperature (20-25°C) before use.
2. **Wash Buffer** - Dilute 30mL of 20x Wash Buffer Concentrate with 570mL of deionized or distilled water to prepare 600mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.
3. **Standard/Sample (2 vials)** - Rabbit IL-12 Standard has a total of 2 vials. Each vial contains the standard sufficient for generating a standard curve. Reconstitute the Standard with 1.0mL of Standard/Sample

Diluent. This reconstitution produces a stock solution of 1600 pg/mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 500µL of Standard/Sample Diluent into 800pg/mL tube and the remaining tubes. Use the stock solution of 1600pg/mL to produce a 2-fold dilution series (below). Mix each tube thoroughly (vortex 20 sec for each of dilution step) and change pipette tips between each transfer. The 1600 pg/mL standard serves as the high standard. The Standard/sample Diluent serves as the zero standard (0 pg/mL).

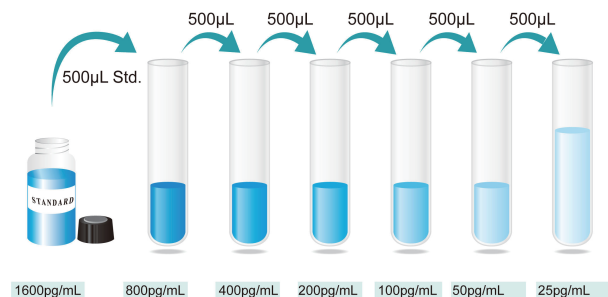

Preparation of Rabbit IL-12 standard dilutions

**\*If you do not run out of re-melting standard, store it at  $-20^{\circ}\text{C}$ . Diluted standard shall not be reused.**

4. **Working solution of Biotin-Conjugate anti-Rabbit IL-12 antibody(1 vials)**-The lyophilized Detection Antibody should be stored at 4°C to  $-20^{\circ}\text{C}$  in a manual defrost freezer for up to 6 months, if not used immediately. Centrifuge for 1 min at 6000 x g to bring down the material prior to open the vial. The vial contains sufficient Detection Antibody for a 96-well plate. Add 110µL of sterile Biotin-Conjugate antibody Diluent to each vial and vortex 30 sec to obtain the stock solution. If the entire 96-well plate is used, take 100µL of detection antibody stock solution into 10 mL of Biotin-Conjugate antibody Diluent to make working dilution of Detection